



JOBST®

JoViPak®

Product Catalog
v2



JOBST® JoViPak®

Garments for Nighttime Lymphedema
Management, and Rehabilitation

Total Body Solutions

Compression where you need it with JOBST® JoViPak®

JOBST® JoViPak® garments offer easy-to-use alternatives to nighttime compression bandaging, as well as rehabilitation for injuries. These foam filled garments provide pressure into the tissue and, with the recommended JoViJacket, provide additional compression for optimal therapeutic results. Below are ready-to-wear garments; however JOBST® JoViPak® offers **custom made garments** for **any part of the body**, in **any size**, and with an array of **option** choices.

JOBST® JoViPak® ready-to-wear and made-to-measure leg and arm garments contain the JOBST® JoViPak® garment and the outer JOBST® JoViJacket. Both components together are designed to provide approximately 20-30mmHg compression.

Head & Neck

- Chin Straps
- Face Masks
- Neck Pads
- Eye Pad

Upper Body

- Arm Sleeve
- JoViLiner
- Rotator Cuff
- Cock Up Splint
- Elbow Wrap

Torso

- Bellisse®
- Breast & Chest Wall Pads

Genital Area


- Genital Pads

Lower Body

- Leg Sleeve
- JoViLiner
- BinkyBoot (for infants & small children)
- Thigh Wrap
- Knee Wrap
- Calf Wrap
- Foot and Lower Leg Wrap
- Ankle Wrap

Pads & PitPaks*

- JOBST® JoViPads
- JOBST® PitPaks

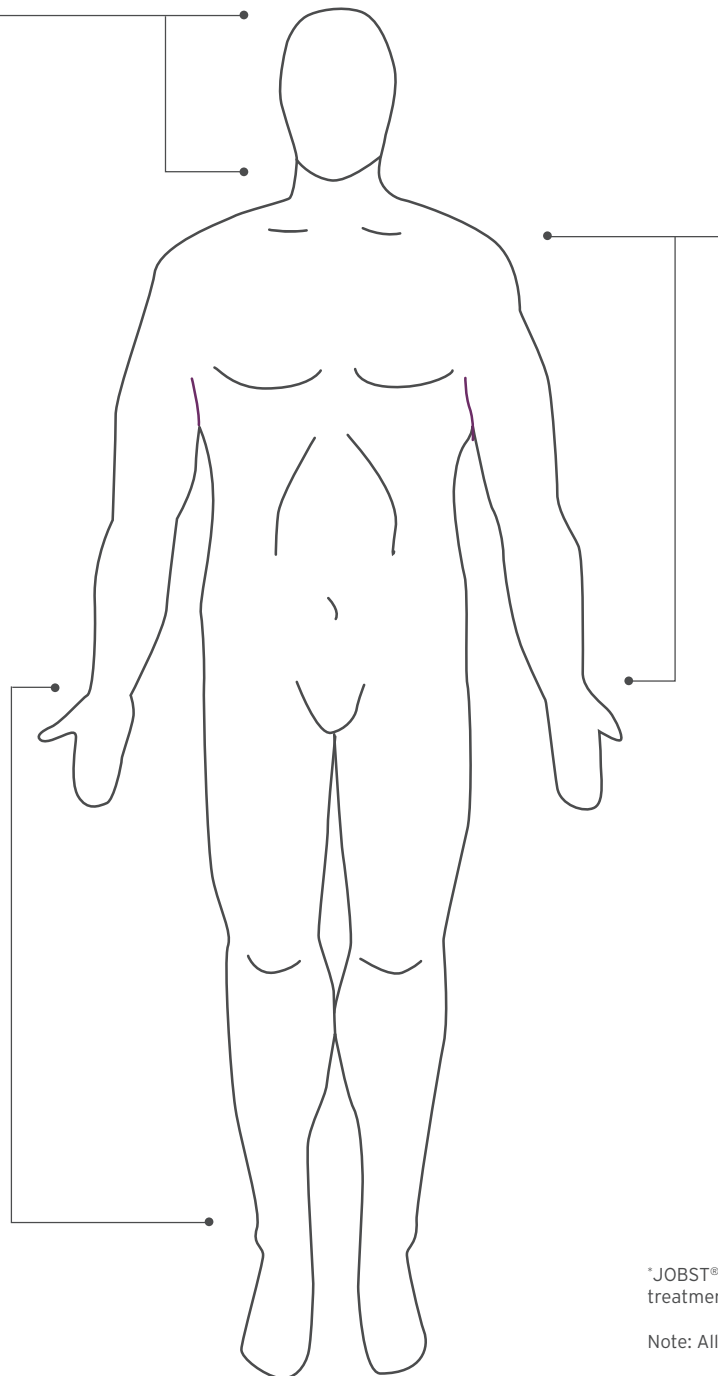

*JOBST® JoViPads and JOBST® PitPaks can be used in treatment on any area of the body.

Note: All product order forms can be found on jovipak.com


Unique Features and Benefits

Continuous Flow channeling technology.

- Directs lymph flow to functional alternate pathways
- Supports lymph fluid drainage

● Polartec® Power Dry® fabric offers special microfiber yarns and a unique weave to disperse moisture and heat across a wide surface area.

- Helps quick drying of fabric

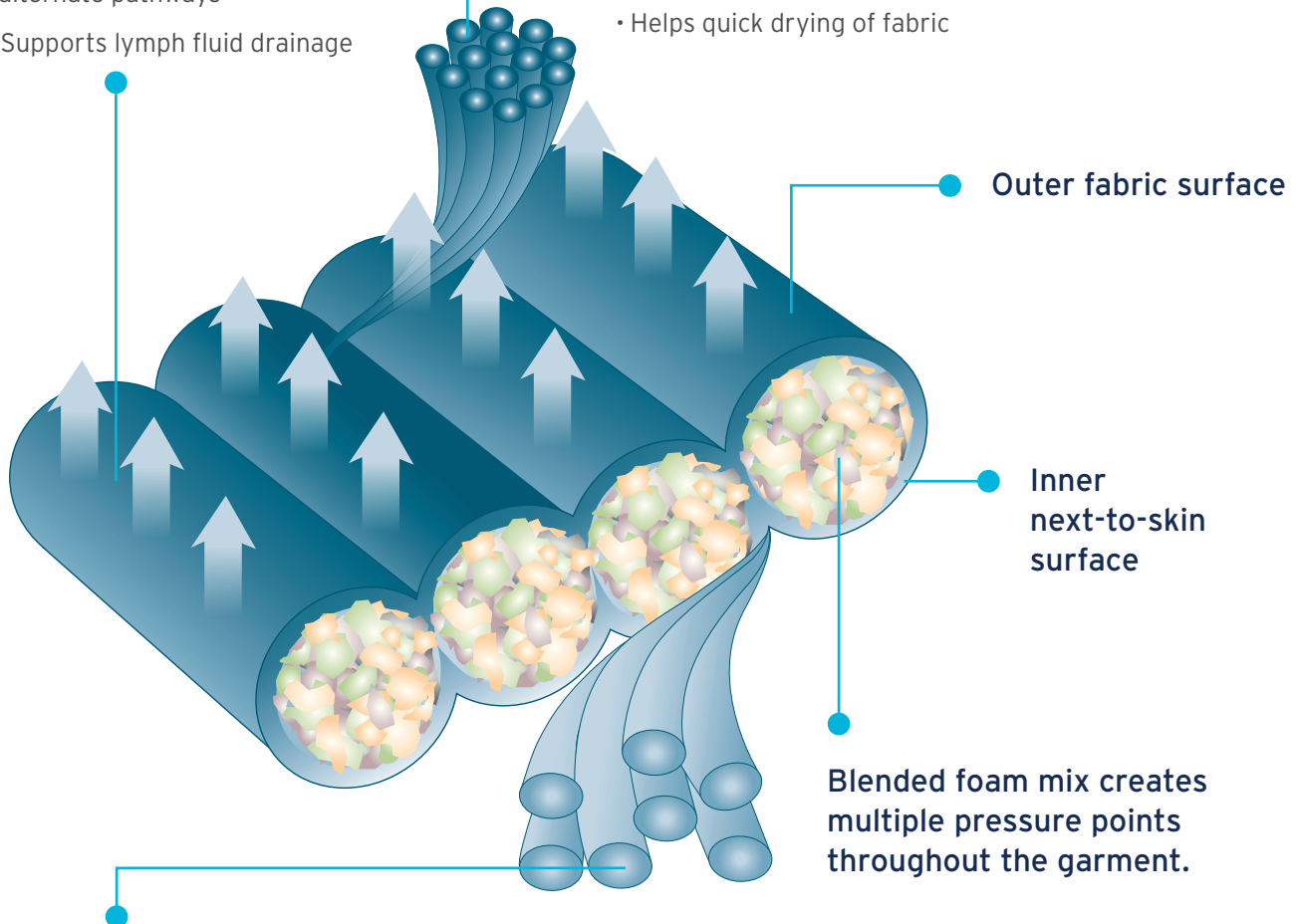

● Outer fabric surface

● Inner next-to-skin surface

● Blended foam mix creates multiple pressure points throughout the garment.

● Smooth, highly durable yarns optimized to help move moisture away from the skin.

- For comfortable wearing

The **Polartec® Power Dry®** fabric is also available in **silkweight**, a thinner version which is ideal for low profile slimline garments. **Silkweight** is often recommended for patients having donning and doffing issues.

Garments made of **Organic Cotton / Spandex** may be a better choice for patients with a skin sensitivity to polyester, and for products with Velcro® closures as Velcro® closures will snag the **Polartec® Power Dry®** fabric.

Table of Contents

READY-TO-WEAR UPPER BODY GARMENTS

Arm

Classic Arm Sleeve AG1	5
Classic Arm Sleeve AG1 JoviJacket	5
JoviLiner Arm	5

Head & Neck

Chin Strap - Standard	7
Chin Strap - Extended	7
Face Mask - Half	7
Face Mask - Full	8
Neck Pad - Universal	8
Neck Pad - Peri-Auricular	8
Head & Neck Extender	9
Eye Pad	9

Torso

Bellisse®	11
Belisse® Extender	11
Drain Pocket	11
Axilla Pad	12
Chest Wall Pocket Pad	12
Cleavage Pad	13
Double Mastectomy Pad	13
Inframammary Pad	13
Lateral Pad	13
Mini-Axilla Pad	14
Padded Insert	14
Post-Lumpectomy Pad	14
Serratus Anterior Pad	15
Unilateral Post-Mastectomy Pad	15

READY-TO-WEAR LOWER BODY GARMENTS

Leg

Classic Lower Leg AD	17
Classic Lower Leg AD JoviJacket	17
JoviLiner Full Leg	17

Genital Area

Female Genital Pad	19
Male Genital Pad	19

JOBST® JOVIPADS, JOBST® PITPAKS & DONNING AIDS

JOBST® JoViPads

Round Pad	20
Oval Pad	20
Dorsum Pad	20
Dorsum Insert	20
Palm Pad	21
Palm Insert	21
Malleolus Pad	21
Kidney Pad	22
Rectangle Pad	22

JOBST® PitPaks

Round PitPak	23
Oval PitPak	23
Crescent PitPak	24
Dorsum PitPak	24
Half Circle PitPak	24
Rectangle PitPak	24
Inframammary PitPak	24
Full Channeled Cylinder PitPak	25
Half Channeled Cylinder PitPak	25
Tapered Strip PitPak	25

Donning Aids

Arion	26
Dycem®	26

CUSTOM-MADE OPTIONS

Upper Body	27
Lower Body	28
Fabrics & Colors	29

REHAB GARMENTS

Rotator Cuff & Rotator Cuff JoViJacket	31
Elbow Wrap	31
Cock-Up Splint	31
Thigh Wrap	32
Knee Wrap - Standard & Extended	32
Calf Wrap	32
Foot & Lower Leg Wrap and Foot & Lower Leg Wrap JoviJacket	33
Ankle Wrap & Ankle Wrap JoviJacket	33



Classic Arm Sleeve AG1

The **Classic Arm Sleeve AG1** was JoViPak's original product and continues to be the most popular.


This gradient pressure arm sleeve is filled with chipped foam. The special foam assists in breaking down trapped proteins to help move lymphatic fluid and help manage lymphatic swelling. The sleeve is engineered to eliminate the need for nighttime bandaging.

The Ready-to-Wear design features a gauntlet style hand and lateral rise at the shoulder. An anatomically positioned thumb stub and an anticubital dart adds comfort.

This product is available with standard channeling in our signature Royal Blue or Black. JOBST® JoViPak® recommends ordering a **Classic Arm Sleeve AG1 JoViJacket** with this garment. JoViJackets provide the additional compression needed for maximum fit & effectiveness.

For alternative styles or colors, please choose a Custom-Made garment.

Fabric	Polartec® Power Dry®	
Color	Royal Blue standard channeling	

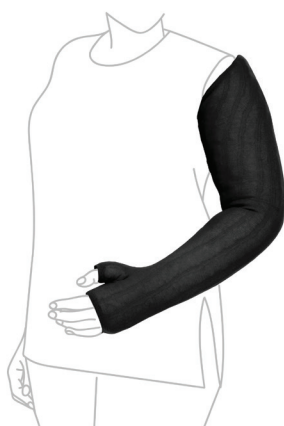

Classic Arm Sleeve AG1 JoViJacket

This outer JoViJacket adds additional compression to the foam filled **Classic Arm Sleeve AG1**. A JoViJacket is recommended with the appropriate garment to enhance the fit and to maximize the therapeutic benefit.

Fabric	Nylon and Spandex	
Color	Black	



JoViLiner Arm

The **JoViLiner Arm** is made of a soft cotton & Spandex blend and is designed to provide comfort and protection for sensitive areas under either short-stretch bandaging or pumps.

Fabric	Organic Cotton and Spandex	
Color	Ivory	


Size Chart

Arm Size	Least Wrist Circumference	Least Elbow Circumference	Axilla Circumference
	"C" Landmark	"E" Landmark	"G" Landmark
Small	15.5 - 16.5	23 - 26	27.5 - 32
Medium	16.5 - 18	26 - 29	32 - 36
Large	18 - 19	29 - 33	36 - 41
XLarge	19 - 20.5	33 - 37.5	41 - 45

Length Chart

Wrist to Axilla (C to G) <small>(measured with arm hanging down at side)</small>	
Regular	43
Long	45.5

Landmark Chart

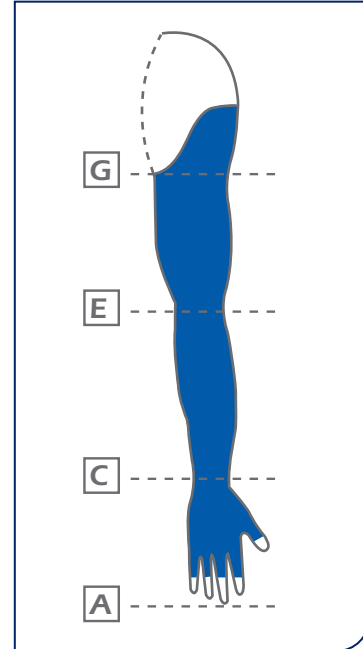

* All measurements are in centimeters (cm)



Chin Strap - Standard

The **Chin Strap - Standard** is designed to address edema of the neck, mandible and the anterior/lateral aspects of the face (partially covers the ears). The open neck seam accommodates a choice of padding. The Velcro® Brand straps at the neck and top of head can be adjusted to patient's comfort level, making the garment suitable for overnight wear.


Fabric	Organic Cotton and Spandex	
Color	Ivory	

Note: this garment does not have any padding. JOBST® JoViPak® recommends a Neck Pad - Universal with this product for maximum effectiveness.



Chin Strap - Extended

Designed to address edema of the neck, mandible and the lateral aspects of the face, the **Chin Strap - Extended** covers the ear, the lateral aspect of the head and behind the ear. Velcro® straps at the neck and top of head (occipital head strap) can be adjusted to patient's comfort level, making the garment suitable for overnight wear.

Fabric	Organic Cotton and Spandex	
Color	Ivory	


Note: this garment does not have any padding. JOBST® JoViPak® recommends a Neck Pad - Peri-Auricular with this product for maximum effectiveness.



Face Mask - Half

The **Face Mask - Half** is a modified version of the Chin Strap. This loose foam filled garment covers the cheek, eye and forehead and has Velcro® straps at the neck and the top and back of head that can be adjusted to patient's comfort level, making the garment suitable for overnight wear.

Fabric	Organic Cotton and Spandex	
Color	Ivory	

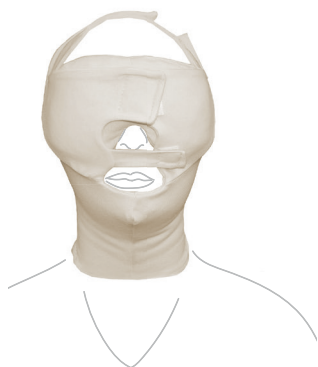

Face Mask - Full

The **Face Mask - Full** is designed to cover the entire face except the nose and lips. This mask has adjustable Velcro® straps on the top and back of the head, at the neck, and above the top lip. The adjustable strap across the upper half of the face allows for full vision when needed.

Fabric	Organic Cotton and Spandex	
Color	Ivory	

Note: this garment does not have any padding. JOBST® JoViPak® recommends a Neck Pad - Universal with this product for maximum effectiveness.



Neck Pad - Universal

The **Neck Pad - Universal** is designed to be used with the **Chin Strap - Standard**, addressing bilateral head & neck edema where there is swelling of the neck in front of the ear(s).

Fabric	Organic Cotton and Spandex	
Color	Ivory	



Neck Pad - Peri-Auricular


Neck Pad - Peri-Auricular addresses bilateral head & neck edema and has extensions for both the anterior and posterior of the ear. This item is recommended with either the **Chin Strap - Extended** or with the **Face Mask - Full**.

Fabric	Organic Cotton and Spandex	
Color	Ivory	

Head & Neck Extender



The **Head & Neck Extender** is medical grade Velcro® Brand fasteners that extends the circumference of the Head & Neck garments and can be used either where the garment fastens at the neck or over the top of the head. This item is very useful if a patient is between sizes of a Chin Strap or Face Mask and/or swelling fluctuates a great deal. The Head and Neck Extender measures 6.5 x 12.7cm.

Fabric	Nylon and Polyester	
Color	Ivory	




Eye Pad

An **Eye Pad** is designed to provide bilateral padding over the orbital area.

Fabric	Polartec® Power Dry®	
Color	Buff	

Landmark Chart

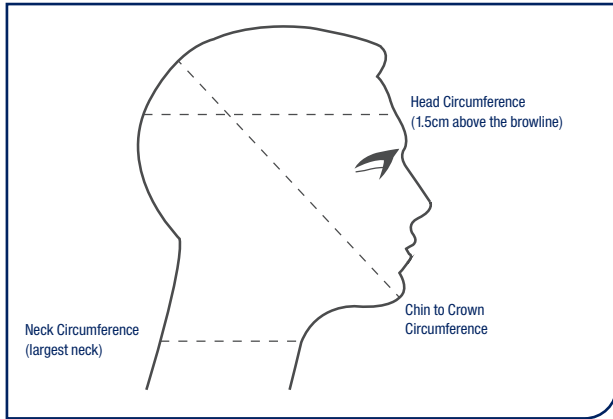

Size Chart

Size	Head Circumference	Neck Circumference	Chin to Crown Circumference
Small	53.3 - 55.9	28 - 33	60 - 63.5
Medium	55.9 - 58.4	33 - 38.1	63.5 - 67.5
Large	58.4 - 61	38.1 - 43.2	67.5 - 71
XLarge	61 - 63.5	43.2 - 48.3	71 - 75
2XLarge (available for Chin Straps only)	63.5 - 66	48.3 - 52.1	75 - 78.5

* All measurements are in centimeters (cm)



Bellisse®

Bellisse® is a medical compression garment that facilitates drainage of excess lymphatic fluid in the chest and breast area following breast cancer and other chest or truncal surgeries.


Fabric	Polyester/Nylon/Spandex			
Color	Black		Buff	



Bellisse® Extender


The **Bellisse® Extender** adds an additional 3 rows adjustability to the rear closure of the Bellisse®. This allows the Bellisse® to be more functional, allow room for padding, and create a more comfortable fit.

Fabric	Organic Cotton and Spandex			
Color	Black		Buff	

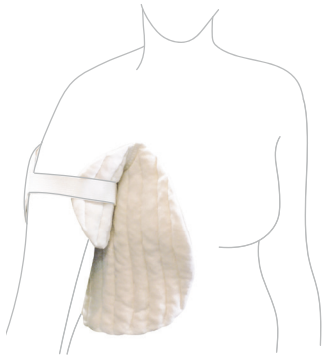


Axilla Pad

The **Axilla Pad** was created in response to a request by a breast surgeon looking for a "tapeless solution" immediately following axillary node dissection. The flat pad is placed against the chest wall and held in place by the Bellisse®.

Fabric	Polartec® Power Dry® with elastic band			
Color	Black		Buff	



Chest Wall Pocket Pad

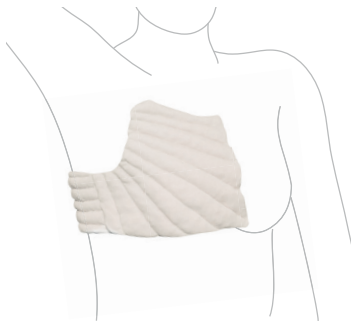


The **Chest Wall Pocket Pad** is a comfortable foam filled pad that conveniently fits into the Bellisse® inner pocket & covers the affected area completely, protecting sensitive tissue. It directs lymph flow around the torso & across the spine, toward the unaffected axilla.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	



Cleavage Pad

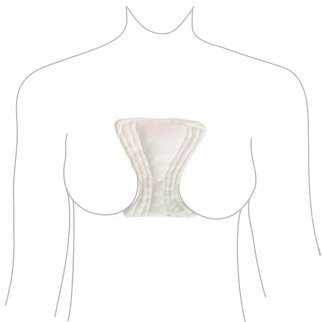

The **Cleavage Pad** is particularly useful under the Bellisse® and mastectomy bras when there is swelling between the breasts. It is also effective for scar management following open heart surgery.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	





Double Mastectomy Pad

The **Double Mastectomy Pad** is a wrap-around torso pad and is secured with Velcro®. It addresses edema on the anterior chest wall as well as the lateral & posterior chest wall areas.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	



Inframammary Pad

The **Inframammary Pad** is an elliptical-shaped pad designed to extend from the serratus anterior to the xiphoid process. The design provides gentle pressure along the inferior breast when used under a Bellisse® or other compression garment.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	



Lateral Pad

The **Lateral Pad** fits into the inside pocket of the Bellisse® and covers the Bellisse® or other compression garment. It is ideal for softening fibrotic areas on the lateral breast.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	





Mini-Axilla Pad

The **Mini-Axilla Pad** has loose fill foam that provides compression into the axilla with a small continuous flow extension. This item should be secured in place by Bellisse® or other compression garment.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	




Padded Insert



Designed to equalize pressure over the site of a mastectomy surgery, a foam filled **Padded Insert** can provide extra comfort and promote a normal silhouette after mastectomy surgery.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	

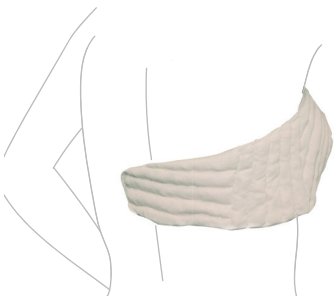

Post-Lumpectomy Pad

The **Post-Lumpectomy Pad** features a posterior extension that addresses swelling and/or radiation fibrosis over the affected breast and serratus anterior. This pad covers the affected breast completely and extends to the spine.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	





Serratus Anterior Pad

Serratus Anterior Pads direct lymph flow laterally down to the inguinals and posteriorly towards the spine. This pad can be used with the Bellisse® or most mastectomy bras.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	




Unilateral Post-Mastectomy Pad

Unilateral Post-Mastectomy Pads are for mastectomy patients with chest wall edema and/or radiation fibrosis over the chest wall and the serratus anterior, with normal lymph drainage on the unaffected side.

Fabric	Polartec® Power Dry®			
Color	Black		Buff	



Landmark Chart

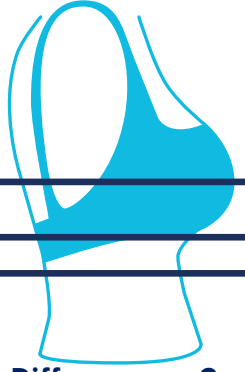

Enter Band Measurements Here

Step 2

Step 1a

Step 1b

Difference	Cup Size
2" or less	A/B
2" - 4"	C/D
4" - 6"	DD/E

If more than 6" go up one band size.

Bellisse Size:



Step 1: Determine the band size:

- Measure circumferentially below the bust where a regular bra sits.
- Drop the tape 2" - 3" and measure around the bottom of the rib cage.
- Use the larger of the two rib cage measurements as the band size. Round odd numbers up.

Step 2: Determine the cup size:

- Measure circumferentially around the fullest part of the bust.
- Subtract the larger ribcage measurement from the bust measurement to determine cup size.

Size Chart - Breast & Chest Wall Pads

Size	Cup Size	Rib Cage Circumference
Small	A/B Cup	32" - 36" (81.3 - 91.4cm)
Medium	C Cup	38" - 40" (96.5 - 101.6cm)
Large	D/DD Cup	42" - 44" (106.7 - 111.8cm)
XLarge	E Cup	46" - 50" (116.8 - 127cm)



Classic Lower Leg AD

The **Classic Lower Leg AD** is for patients who have swelling in their leg and foot as a result of lymphedema or chronic venous insufficiency. The sole of the foot extends to the tip of the longest toe, with the top cut back at an angle to expose the toes. The Ready-to-Wear design is a cost-effective tool to assist in the management phase of lymphedema.

JOBST® JoViPak® recommends ordering the **Classic Lower Leg AD JoViJacket** with this garment to provide the additional compression needed for maximum fit & effectiveness.

Fabric	Polartec® Power Dry®			
Color	Black		Royal Blue	



Classic Lower Leg AD JoViJacket

A **Classic Lower Leg AD JoViJacket** provides the compression needed for nighttime lymphedema management and has a non-slip sole to offer extra safety.

A JoViJacket is recommended with the appropriate garment to enhance the fit and to maximize the therapeutic benefit.


Fabric	Nylon and Spandex	
Color	Black	



JoviLiner Full Leg

The **JoViLiner Full Leg** has continuous flow channels from the tip of the toes to the groin and covers to the base of the toenail beds with stitching that divides all but the third and fourth toes. It is intended for use under short-stretch bandaging or with pumps.

JoViLiners are made of a soft cotton & Spandex blend, designed to provide comfort and protection to sensitive areas.

Fabric	Organic Cotton and Spandex	
Color	Ivory	


Size Chart

Leg Size	Least Ankle Circumference	Widest Calf Circumference	Least Knee Circumference
	"B" Landmark	"C" Landmark	"D" Landmark
Medium	22 - 23	38 - 40.5	34.5 - 37
Large	23 - 24	40.5 - 43	37 - 39.5
XLarge	24 - 25.5	43 - 45.5	39.5 - 42
2XLarge	25.5 - 27	45.5 - 48	42 - 44.5

Length Chart

Floor to Tibial Tuberosity (A to D) (measured medially)	
Regular	42

Landmark Chart

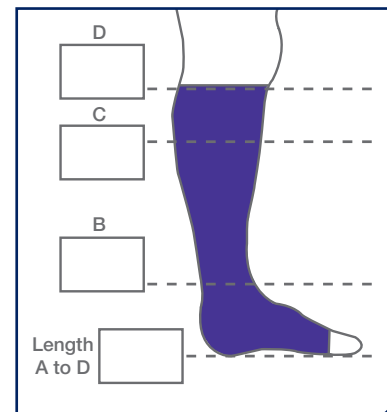


* All measurements are in centimeters (cm)

Female Genital Pad



The bikini cut of the **Female Genital Pad** is designed so that the anterior width covers most hysterectomy scars. This item is a foam filled liner with a vertical continuous flow stitching pattern that directs the flow of lymph fluid away from the pubis and toward axillary collection areas. Use this product under a JOBST® compression garment for protection, added comfort and optimal results.

Shown here is the **Female Genital Pad** with straps (containing latex) over the iliac crest to hold the garment in place. This product is also available without straps, but would always need to be held in place by another garment.


Fabric	Polartec® Power Dry®	
Color	Buff	
Sizing	Available in 5 sizes determined by pant size: X-Small (0-4), Small (4-8), Medium (10-12), Large (14-16), and X-Large (18-22)	

Male Genital Pad




The **Male Genital Pad** is recommended for scrotal swelling and/or fibrosis. The **Male Genital Pad** is a foam filled liner with a vertical continuous flow stitching pattern that directs the flow of lymph fluid away from the testicular area and toward an alternate functional lymphatic pathway. Use this product under a JOBST® compression garment for protection, added comfort, and optimal results.


Shown here is the **Male Genital Pad** with straps (containing latex) over the iliac crest to hold the garment in place. This product is also available without straps but would always need to be held in place by another garment.

Fabric	Polartec® Power Dry®	
Color	Buff	
Sizing	Available in 5 sizes determined by scrotal circumference: X-Small (8"), Small (11"), Medium (13"), Large (16.5"), and X-Large (20.5")	

Round Pad

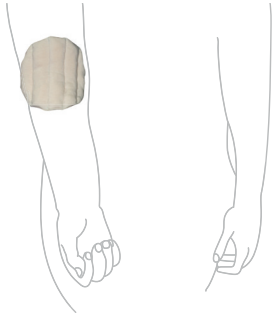


This **Round Pad** provides protection from bandaging in the delicate antecubital fossa area but can also be used in any small area where fibrosis is a problem. Its size allows it to be easily used under bandages, pump appliances or compression garments.


Fabric	Organic Cotton and Spandex	
Color	Ivory	
Sizing	Available in 1 size: 12.5 cm	



Oval Pad

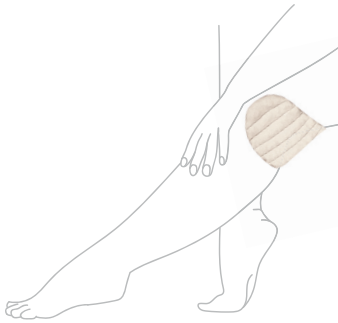


Oval Pads provide protection from bandaging in the delicate popliteal fossa area. The **Oval Pad** can be used in any small area where fibrosis is a problem and fits easily under bandaging, pump appliances or compression garments.


Fabric	Organic Cotton and Spandex	
Color	Ivory	
Sizing	Available in 1 size: 12.5 x 20 cm	



Dorsum Pad

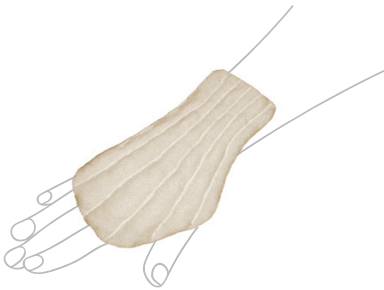

The **Dorsum Pad** is designed to fit over the dorsum of the hand. Stitched channels follow the metacarpals, with finely chipped foam sinking down into the spaces between. It is often combined with the **Palm Pad** under bandaging, whenever hand swelling is a problem.


Fabric	Organic Cotton and Spandex	
Color	Ivory	
Sizing	Available in 1 size: 5.5 x 4.75 cm	





Palm Pad

The **Palm Pad** is designed to help maintain proper functional arches of the hand under bandage pressure.


Fabric	Organic Cotton and Spandex	
Color	Ivory	
Sizing	Available in 1 size: 8.58 x 10 cm	



Malleolus Pad



Medial and/or lateral anatomical indentations around the malleoli can be filled by the **Malleolus Pad**. Using this item under bandaging or a JOBST® JoViPak® nighttime compression garment can provide gentle compression to an area that may not otherwise be addressed.


Fabric	Organic Cotton and Spandex	
Color	Ivory	
Sizing	Available in 1 size: 17.5 cm	





Kidney Pad


As with all JOBST® JoViPads, **Kidney Pads** can be placed in a variety of areas. Placing this item on the lateral chest at the axilla, on the distal anterior of the knee, or on the groin at the inner thigh are examples for this versatile product.


Fabric	Organic Cotton and Spandex	
Color	Ivory	
Sizing	Available in 3 sizes: 7.5 x 17.5 cm, 12.5 x 25 cm, 17.5 x 32.5 cm	



Rectangle Pad

Rectangle Pads are highly versatile, designed to be worn almost anywhere additional pressure is needed to provide better therapeutic results. Consider this product for the lateral chest wall, anterior tibialis, lateral aspect of the forearm, etc.

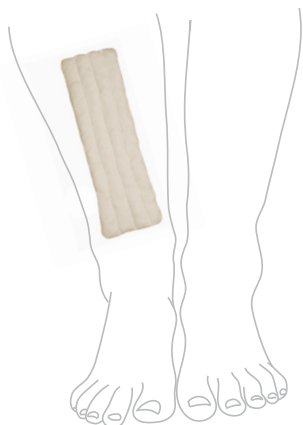


Fabric	Organic Cotton and Spandex	
Color	Ivory	
Sizing	Available in 2 sizes: 10 x 20 cm, 10 x 30 cm	



Round PitPak*




Round PitPaks have multiple uses where progressive tissue change is needed. The small size of this product allows it to be used in smaller areas where there is a fibrotic scarring or radiated tissue.

Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in Loose fill 10 cm or 3 Channeled 15 cm	

Oval PitPak*




A versatile and simple item to use, the **Oval PitPak** could be considered for multiple places that need special attention (i.e., the ulnar aspect of the forearm, etc.),

Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in Loose fill 10 x 15 cm or 3 Channeled 13.9 x 21.5 cm	

Crescent PitPak*



Crescent PitPaks follow the natural curves of the body and address stubborn woody fibrosis in such areas as the base of the toes and the inferior aspect of the breast.

Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in Loose fill 3.8 x 10 cm or 3 Channeled 5 x 12.5 cm	




*Note: JOBST® PitPaks can be part of an aggressive course of therapy and should only be undertaken with the supervision of a qualified therapist or similar medical professional.

Dorsum PitPak*




The **Dorsum PitPak** should be used with a short stretch bandage to soften long standing or woody fibrosis on the dorsum of the hand.

Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in 1 size: 10 x 12.5 cm	

Half Circle PitPak*



A **Half Circle PitPak** consists of multiple rows of cherry pit channeling and is often ideal for use in the sub-mandibular region on patients with head & neck edema.


Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in 1 size: 15 x 20 cm	



Rectangle PitPak*



The **Rectangle PitPak** features multiple rows of cherry pit channeling and may be useful in a variety of areas, providing massaging action as the patient moves.


Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in 1 size: 15 x 22.5 cm	



Inframammary PitPak*




The **Inframammary PitPak** is an elliptical-shaped pad of single row cherry pits for the inferior aspect of the breast.

Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in 4 sizes: 10 x 21.5 cm, 11.43 x 22.5 cm, 12.5 x 24.1 cm, 14 x 25 cm	




*Note: JOBST® PitPaks can be part of an aggressive course of therapy and should only be undertaken with the supervision of a qualified therapist or similar medical professional.

Full Channeled Cylinder PitPak*

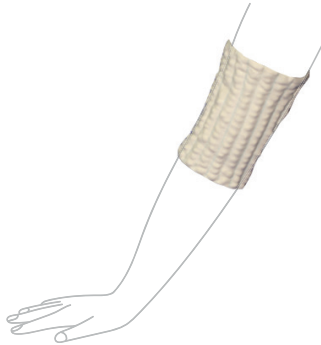

Full Channeled Cylinder PitPaks are used under short-stretch bandaging to soften circumferential woody or long-standing fibrosis. This product has multiple single rows of cherry pits sewn around the entire garment.


Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in 4 sizes: 17.5-22.5 x 20 cm, 22.5-30 x 20 cm, 27.5-37.5 x 28 cm, 35-42.5 x 28 cm	



Half Channeled Cylinder PitPak*



Half Channeled Cylinder PitPaks are used under short-stretch bandaging to soften circumferential woody or long-standing fibrosis. Multiple single rows of cherry pits are sewn into half of the width of this product.


Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in 2 sizes: 18 x 25.5 cm, 20 x 30.5 cm	



Tapered Strip PitPak*



The **Tapered Strip PitPak** contains two channels of loose fill cherry pits comprise this garment which is often ideal for the lateral forearm.

Fabric	Organic Cotton and Spandex with 100% cleaned & sterilized cherry pits	
Color	Ivory	
Sizing	Available in 1 size: 5.71 x 17.5 cm	




*Note: JOBST® PitPaks can be part of an aggressive course of therapy and should only be undertaken with the supervision of a qualified therapist or similar medical professional.



Arion

Arion donning aids are made from very smooth material, with a special coating which together with double-layer product design create an extremely smooth gliding product, this minimizes the friction between the compression garment and skin. Arion Easy-Slide is user-friendly for open toe compression garments.

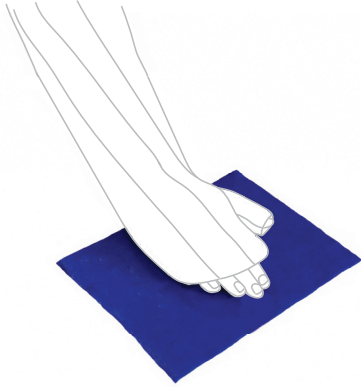

Dycem®

Dycem® facilitates easier donning of both upper and lower extremity garments. Dycem is placed on the floor or a table top to provide a non-slip surface that the patient can push against as they don the garment. Dycem® can be easily cleaned with soap and water to renew its non-slip surface.

JOBST® JoViPak® - A Solution For Every Need

JOBST® JoViPak® offers Custom-Made garments for nighttime compression to address any part of the body. The wide range of Custom-Made garments include many individual style, fabric, and color options.

Custom-Made Garments - Upper Body		
Arm Sleeve AG JoViJacket	Arm Sleeve AC1	InnaSleeve (with SUPER Powernet)
Arm Sleeve AG1	Arm Sleeve AC1 JoViJacket	Shoulder Torso Arm Sleeve
Arm Sleeve AG1 JoViJacket	Arm Sleeve CG1	Shoulder Torso Arm Sleeve JoViJacket
Arm Sleeve AD	Arm Sleeve CG1 JoViJacket	Vest - Half Padding
Arm Sleeve AD JoViJacket	Busti	Vest JoViJacket
Arm Sleeve AE	Busti JoViJacket	



Arm Sleeve AG1



Arm Sleeve AG1 with recommended JoViJacket



Busti



InnaSleeve



Shoulder Torso Arm Sleeve (with optional Padded Torso and One Piece Arm Sleeve)



Vest

Custom-Made Garments - Upper Body		
Arm Sling	Padded Torso - Vertical	Two Piece Arm Sleeve
Crotch Strap (option for Custom Made Vest only)	Padding To Hemline	Two Piece Arm Sleeve JoViJacket
Dorsum Pad	Palm Pad	Zipper - Dorsum to Forearm
Full Padding	Prepaid Reduction	Zipper - Elbow to Axilla
Padded Torso - Horizontal	Stitched Finger Glove	Zipper - Wrist to Elbow



Two Piece Arm Sleeve



Shoulder Torso Arm Sleeve with optional Bilateral Arms, Padded Torso, Arm Sling and Dorsum Zipper



Arm Sleeve AG with optional Dorsum to Mid-Forearm Zipper, Elbow to Axilla Zipper and Stitched Finger Glove

Custom-Made Garments - Lower Body		
ADVI	Hip Hugger Full Leg	Leg AD JoViJacket
ADVI JoViJacket	InnaBoot Lower Leg (with SUPER Powernet)	Leg AF1
Boxer	InnaBoot Full Leg (with SUPER Powernet)	Leg AF1 JoViJacket
Boxer JoViJacket (SUPER Powernet)	Leg AG	Leg AB1
Boxer Capri	Leg AG JoViJacket	Leg AB1 JoViJacket
Boxer Capri JoViJacket (SUPER Powernet)	Leg AG1	Leg DG1
Hip Hugger	Leg AG1 JoViJacket	Leg DG1 JoViJacket
	Leg AD	



ADVI



Boxer with recommended JoViJacket



Boxer Capri



Hip Hugger (posterior)



InnaBoot Lower Leg



Leg AD



Leg AG1

Custom-Made Garments - Lower Body		
ADVI Foot Style	Malleolus Pad - Medial	Zipper - Ankle to Knee
Donning Loops	Pull Tabs (option for Custom-Made InnaBoots only)	Zipper - Knee to Groin
Malleolus Pad - Lateral	Reduction Option	



ADVI Foot Style



Pull Tabs (InnaBoot only)



Zipper - Knee to Groin

Available Custom-Made Fabrics & Colors

Upper Body

CUSTOM ARMS		
Organic Cotton and Spandex <small>(InnaSleeve in Cotton and Spandex only)</small>		
● Black	● Ivory	● Royal Blue
Polartec® Power Dry®		
● Black	● Buff	● Navy Blue
● Pink	● Royal Blue	○ White <small>(soft pink hue)</small>
● Stainless Steel	● Plum	
Polartec® Silkweight		
● Black	● Blue Ridge	
JOVIJACKETS*		
● Black	○ White	

CUSTOM SHOULDER-TORSO ARM SLEEVES		
Polartec® Power Dry®		
● Black	● Buff	○ White <small>(soft pink hue)</small>
● Pink	● Royal Blue	● Navy Blue
● Stainless Steel	● Plum	
Polartec® Silkweight		
● Black	● Blue Ridge	
JOVIJACKETS*		
● Black	○ White	

CUSTOM BUSTI	
Polartec® Power Dry®	
● Black	● Buff
● Pink	● Royal Blue
● Stainless Steel	● Navy Blue
○ White <small>(soft pink hue)</small>	● Plum
JOVIJACKETS*	
● Black	○ White

CUSTOM VEST		
Organic Cotton and Spandex		
● Black	● Ivory	● Royal Blue
JOVIJACKETS*		
● Black	○ White	

Lower Body

CUSTOM LEGS		
Organic Cotton and Spandex <small>(InnaBoot in Cotton and Spandex only)</small>		
● Black	● Ivory	● Royal Blue
Polartec® Power Dry®		
● Black	● Buff	● Navy Blue
● Pink	● Stainless Steel	● Royal Blue
○ White <small>(soft pink hue)</small>	● Plum	
JOVIJACKETS*		
● Black	○ White	

CUSTOM ADVIS	
Organic Cotton and Spandex	
● Black	● Ivory
CUSTOM ADVI JOVIJACKETS	
● Black	○ White

CUSTOM HIP HUGGERS	
Organic Cotton and Spandex	
● Black	● Ivory
● Royal Blue	
CUSTOM HIP HUGGER JOVIJACKETS*	
● Black	○ White

CUSTOM BOXERS & BOXER CAPRIS	
Polartec® Power Dry®	
● Black	○ White <small>(soft pink hue)</small>
● Stainless Steel	● Buff
● Plum	● Navy Blue
● Pink	● Royal Blue
JOVIJACKETS* (SUPER POWERNET)	
● Black	● Buff
○ White	

* JoViJackets are external compression garments recommended for use with matching products to maximize therapeutic benefit.




JOBST® JOVIPAK® Rehab Garments for a Faster Recovery

Whether getting athletes back in the game quicker or knee replacement patients back on their feet sooner, the JOBST® JoViPak® line of rehabilitation products can be an invaluable tool.

- Reduce swelling associated with inflammation
- Provide compression and soft splinting
- Support and protect joints
- Promote circulation and gently massage soft tissues

JOBST® JoViPak® Rehab Garments feature:

- Highly resilient shredded foam which creates high and low pressure points in the skin and underlying tissue
- Improved micro-circulation through gentle compression
- Adjustable and washable material, may help improve patient compliance through easy and independent use

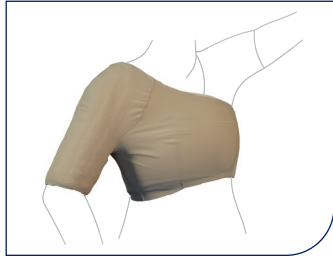

Rotator Cuff and Rotator Cuff JoViJacket

The **Rotator Cuff** garment is designed to provide comfort and support for the rotator cuff. This product is extremely helpful in stabilizing AC separation (shown with the optional, recommended, JoViJacket).

Rotator Cuff

Size Chart*

Size	Circumference (cm)	
	Widest Bicep	Chest
Small	25.5 - 30.5	86.5 - 94
Medium	30.5 - 35.5	97 - 103
Large	35.5 - 40.5	103 - 112
XLarge	40.5 - 45.5	112 - 124


Colors	
Black	
Buff	

Rotator Cuff JoViJacket

Size Chart*

Size	Circumference (cm)	
	Widest Bicep	Chest
Small	25.5 - 30.5	86.5 - 94
Medium	30.5 - 35.5	97 - 103
Large	35.5 - 40.5	103 - 112
XLarge	40.5 - 45.5	112 - 124

Colors	
Black	
Buff	

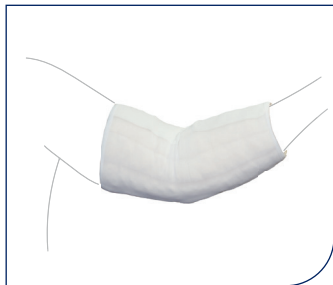

Elbow Wrap

The **Elbow Wrap** is recommended for edema and pain management in the elbow. This wrap allows easy adjustability and functional elbow flexion to approximately 90°. The garment features an attached Nylon and Spandex outer wrap for effectiveness.

Size Chart*

Size	Elbow Circumference (cm)
Small	25.5 - 30.5
Medium	30.5 - 35.5
Large	35.5 - 40.5
XLarge	40.5 - 45.5

Colors	
Black	
Ivory	



Cock-Up Splint


The easy design of the **Cock-Up Splint** combines the features of both a hard and soft splint. It wraps at the wrist and the forearm and incorporates a removable metal splint.

Size Chart*

Size	Widest Palm @ Thumb Web Circumference (cm)
Small	16.5 - 19
Medium	19 - 21.5
Large	21.5 - 24
XLarge	24 - 26.5

Colors	
Black	
Ivory	



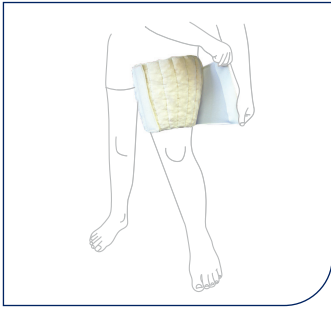


Thigh Wrap

The **Thigh Wrap** extends from lower thigh to groin with adjustable Velcro® Brand closure. The attached Nylon and Spandex wrap-around Powernet® eliminates the need for short-stretch bandaging.

Size Chart*

Size	Mid-Thigh Circumference (cm)
Small	41 - 49.5
Medium	49.5 - 58.5
Large	58.5 - 67.5
XLarge	67.5 - 76

Colors	
Black	■
Ivory	■

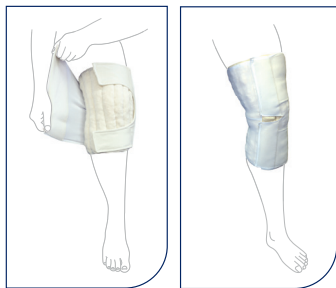

Standard

Extended

Knee Wrap - Standard and Extended

This wrap-style knee garment covers the leg from below the mid-calf to the lower thigh. The Knee Wrap Extended is also available and covers the leg from below the mid-calf to the upper thigh. Either Knee Wrap is ideal for acute injuries or knee replacements. They are designed for acute or traumatic edema when either donning is a concern or where wound dressing is a consideration.

Knee Wrap - Standard

Size Chart*

Size	Knee Circumference (cm)
Small	29 - 34
Medium	34 - 39.5
Large	39.5 - 44.5
XLarge	44.5 - 51


Colors	
Black	■
Ivory	■

Knee Wrap - Extended

Size Chart*

Size	Circumference (cm)	
	Knee	Mid-Thigh
Small	29 - 34	41 - 49.5
Medium	43 - 39.5	49.5 - 58.5
Large	39.5 - 44.5	58.5 - 67.5
XLarge	44.5 - 51	67.5 - 76

Colors	
Black	■
Ivory	■

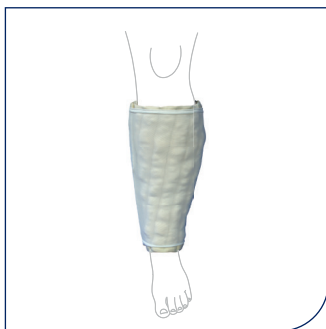

Calf Wrap

The **Calf Wrap** is an adjustable foam-filled wrap-around garment that covers the calf from ankle to knee. The garment is easy to don and doff and features an attached Nylon and Spandex outer wrap for effectiveness.

Size Chart*

Size	Widest Calf Circumference (cm)
Small	33 - 38
Medium	38 - 43
Large	43 - 48
XLarge	48 - 53

Colors	
Black	■
Ivory	■



Foot & Lower Leg Wrap and Foot & Lower Leg Wrap JoviJacket

The **Foot and Lower Leg Wrap** has a unique adjustable design for acute injuries allows for easy donning & adjusts to patient's tolerance. The garment fits like a boot extending to the base of the knee with a full-length Velcro® Brand closures that can be tightened as swelling is reduced. It comes with an inner anterior tibialis pad which provides additional protection over the ankle and settles into the anatomical indentations around the malleoli.

Foot & Lower Leg Wrap Size Chart*

Size	Circumference (cm)	
	Ankle	Widest Calf
Small	19 - 21.5	33 - 38
Medium	21.5 - 24	38 - 43
Large	24 - 26.5	43 - 48
XLarge	26.5 - 29	48 - 53

Colors	
Black	■
Ivory	□

Foot & Lower Leg Wrap JoViJacket Size Chart*

Size	Circumference (cm)	
	Ankle	Widest Calf
Small	19 - 21.5	33 - 38
Medium	21.5 - 24	38 - 43
Large	24 - 26.5	43 - 48
XLarge	26.5 - 29	48 - 53

Colors	
Black	■
White	□



Ankle Wrap and Ankle Wrap JoViJacket

The **Ankle Wrap's** unique adjustable design allows for easy donning and is adjustable per the patient's tolerance. In a rehab setting, the garment can be used to help heal acute injuries (e.g., acute sprains, fractures, and post surgical management).

Ankle Wrap Size Chart*

Size	Ankle Circumference (cm)
Small	19 - 21.5
Medium	21.5 - 24
Large	24 - 26.5
XLarge	26.5 - 29

Colors	
Black	■
Ivory	□

Ankle Wrap JoViJacket Size Chart*

Size	Ankle Circumference (cm)
Small	19 - 21.5
Medium	21.5 - 24
Large	24 - 26.5
XLarge	26.5 - 29

Colors	
Black	■
White	□

Rehab garments are made of Organic Cotton and Spandex blend. Rehab garment JoviJackets are Nylon and Spandex (standard) Powernet, except for the Rotator Cuff JoViJackets which are Nylon, Spandex and SUPER Powernet.

NOTES

JOBST® Educational Resources

Schedule a Virtual In-Service to learn all about JOBST® products! Whether you are new to JOBST® or a longtime friend, we are dedicated to educating about lymphedema and how our products are a crucial part of decongestion and maintenance. Earn 1 Continuing Education Unit (CEU) for 1 hour of participation.



All you need is a computer with audio and a desire to provide the best products to help your patient's manage their lymphedema. With your input our technical advisors will customize the presentation to your needs and your schedule.

JOBST® Compression Institute

The JOBST® Compression Institute serves as a resource for healthcare providers and clinicians seeking further education about vein diseases and lymphedema. This website provides **order forms, literature, videos,** and **links** to our national affiliates. Additionally, the JOBST® Compression Institute offers a variety of nationwide workshops and continuing education classes to train you on JOBST® products and proper measuring and fitting techniques. Up to 9 hours of CEU credits available. Register at www.jobstcompressioninstitute.com.





VELCRO® is a registered trademark of Velcro BVBA.
Please refer to the product label and/or package insert for full instructions on the safe use of these products.



JOBST®,
an Essity brand



/JOBSTUSA



@JOBST_USA



@JOBSTforUSA



jobst-usa.com

BSN Medical Inc., an Essity company
5825 Carnegie Blvd., Charlotte, NC 28209-4633
Tel. (+1) 866 888 5684 Fax (+1) 877 760 4943

63646 R1 ©2020 BSN Medical Inc. G20